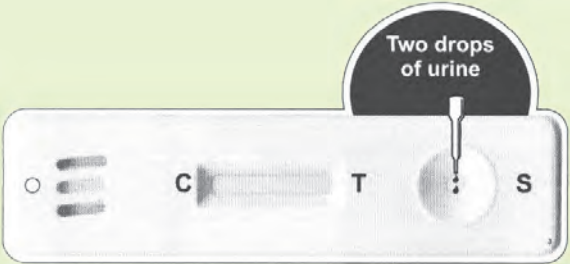

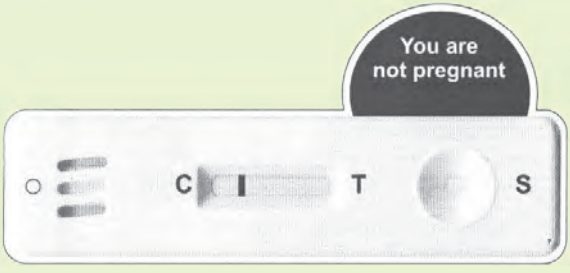



Annexe 2: Instructions for Pregnancy Test using Nischay Kit

The Nischay Kit contains the following:

1. A test card
2. A disposable dropper
3. A moisture absorption packet (not required for testing)

	<ul style="list-style-type: none"> • Collect the morning urine in a clean and dry glass or in a plastic bottle. • Take two drops of urine in the sample well. • Wait for 5 minutes.
	<ul style="list-style-type: none"> • If two violet colour lines come in the test region (T), the woman is pregnant. • If she wants to continue with the pregnancy, advise her to undergo antenatal care. • If she does not want to continue with the pregnancy this time, advise her for safe abortion.
	<ul style="list-style-type: none"> • If the violet colour line in the test region (T) is one only, the woman is not pregnant. • Tell her about family planning methods and help her in choosing the most appropriate one.
	<p>If there is no colour line in the test region (T), repeat the test next morning using a new Pregnancy Test Card.</p>